

MEDICAL

HISTOACRYL[®]

Topical Skin Adhesive



WHY HISTOACRYL[®] FLEXIBLE?

- Histoacryl Flexible has been designed to combine flexibility, strength, and fast polymerization times, to close lacerations and surgical incisions.
- Tensile strength is significantly higher for Histoacryl Flexible, than the leading competitor.
- Close up to 25 linear centimeters with one applicator of adhesive.
- Only one, thin layer of adhesive is needed for full strength closure.
- Flexibility of polymerized adhesive is critical to good cosmetic results.
- Beveled applicator tip provides enhanced control and improved application, for longer, wider incisions.

HISTOACRYL COMPARISON

	Histoacryl Blue & Transparent	Histoacryl Flexible
Composition	n-butyl 2 -cyanoacrylate	n-butyl 2 -cyanoacrylate w/ Plasticizers
Color	Blue & Clear	Blue
Indications	To hold closed clean surgical incisions and cleansed trauma induced lacerations	To hold closed clean surgical incisions and cleansed trauma included lacerations
Color Code <small>(Applicator Packaging)</small>	White and Orange	White and Black
Storage	24 months @ 72°F	24 months @ 77°F
NSN#	6510-01-586-6762 (Histoacryl Transparent) 6510-01-585-2101 (Histoacryl Blue)	N/A



APPLICATOR TIP
Only available with Histoacryl Flexible

1051260-P
Histoacryl Flexible

TS1050071-FP
Histoacryl Transparent

TS1050044-FP
Histoacryl Blue

To learn more visit www.hbfuller.com/medical

HISTOACRYL® BLUE & TRANSPARENT

- **Applicator**- 10 / box

HISTOACRYL® FLEXIBLE

- **Applicator** - 10 / box
- **Applicator Tips** - 10 / box

- Flexibility of an Octyl, tensile strength of a Butyl
- More than 1200 publications site Histoacryl®
- One Histoacryl® Flexible vial, closes 25 linear centimeters
- Histoacryl® Flexible provides a microbial barrier for 7 days, with intact adhesive film
- Precision applicator with optional beveled tip, for longer, wider applications

MICROBIAL BARRIER

In vitro studies have shown that Histoacryl® Flexible acts as a barrier to microbial penetration as long as the adhesive film remains intact.*


Part #	Description	UOM	#/UOM
1051260-P	Histoacryl Flexible 10 / 0.5 mL applicator 10 Applicator Tips	Box	10
TS1050044-FP	Histoacryl Classic - Blue 10 / 0.5 mL applicator	Box	10
TS1050071-FP	Histoacryl Classic - Transparent 10 / 0.5 mL applicator	Box	10

OPEN



Open the package and take out the Histoacryl Tip.

PREPARE APPLICATOR TIP**

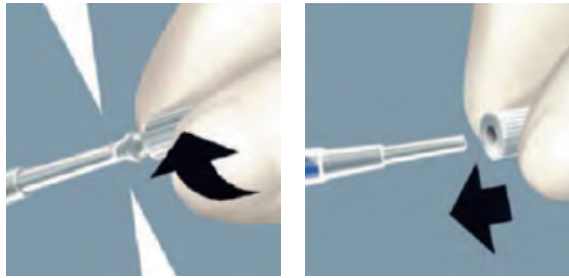

Open the ampoule by twisting off the ribbed tip. Attach the Histoacryl Tip to opened applicator

APPLY



Apply the glue in a thin layer to the approximated wound edges and hold in position for 60 seconds.

The polymerization starts immediately upon contact with the tissue or fluid.

*Histoacryl Flexible inhibits the growth of Staphylococcus aureus, Staphylococcus epidermidis, Escherichia coli, Pseudomonas aeruginosa, Enterococcus faecium, Brevundimonas dimuta and Candida albicans for up to 7 days on inoculated agar plates.

**Applicator Tip only comes with Histoacryl Flex

H.B. Fuller Medical Adhesive Technologies, LLC is the exclusive U.S. partner of Histoacryl Topical Skin Adhesive and B. Braun surgical sutures.